



InSite® | *Listeria* | Rapid Environmental *Listeria* Species Test

Hygiena® InSite® *Listeria* is a quick and easy test for the detection of *Listeria* species in food processing environments. The test contains a liquid media formulated with growth enhancers, selective agents and chromogenic compounds specific to the growth of *Listeria* species.



Benefits:

- Easy to use with no preparation steps required
- Large sponge tip maximizes sample collection for accurate results
- Eliminates costly lab fees for an equivalent test
- Color results are easy to interpret – grey or black indicates presence of *Listeria* species
- Presumptive positive results in as little as 24 hours
- Negative results in 48 hours
- Self-contained tube eliminates risk of contamination


Sensitivity:

- <10 CFU *Listeria* spp.
- Detects heat injured *L. monocytogenes* at low levels (10 – 50 CFU of broth), even in the presence of high numbers of competing organisms

NATURAL LIGHT			
NEGATIVE	<i>LISTERIA</i> SPP.	<i>LISTERIA</i> SPP.	<i>LISTERIA</i> SPP.
—	+ spp.	+ spp.	+ spp.
			
YELLOW	SPOTS	GREY	BLACK
PRESUMPTIVE POSITIVE FOR <i>LISTERIA</i> SPP.			

Product No. (REF)	Description	Quantity
ILC050	InSite <i>Listeria</i> (Environmental <i>Listeria</i> spp.)	50
ILC100	InSite <i>Listeria</i> (Environmental <i>Listeria</i> spp.)	100

InSite® *Listeria* Procedure

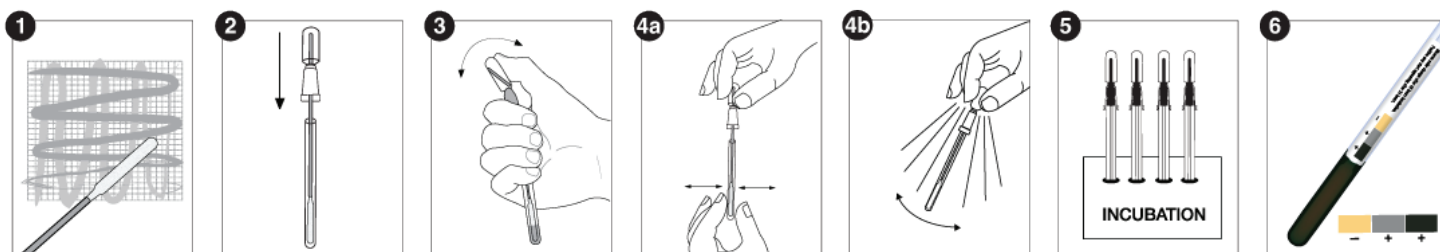

Related Products

Digital Dry Block Incubators

Digital incubators are capable of up to 105 °C and compatible with any Hygiena test device. Available in two sizes. Blocks are sold separately.



InSite® *Salmonella*

InSite *Salmonella* is a quick and easy test for the detection of *Salmonella* species in 24 – 48 hours. The test contains a liquid media formulated with selective agents, growth enhancers, and color changing compounds specific to the growth of *Salmonella* species.



MicroSnap® EB

MicroSnap EB a rapid test for detection and enumeration of *Enterobacteriaceae* (EB) bacteria. Results are available within the same day.



MicroSnap® Coliform & *E. coli*

MicroSnap Coliform & *E. coli* are rapid tests for detection and enumeration of Coliform bacteria. Single figure organisms can be detected in the same day.

MicroSnap® Total

MicroSnap Total delivers total viable count (TVC) results in the same day.

BAX® System Instruments

The BAX System harnesses the power of DNA-based foodborne pathogen detection to help you make product decisions with confidence. An extensive assay portfolio, multiple platforms, and flexible protocols enable food labs of any size to take advantage of this benchmark technology.



Product No. (REF)	Description	Quantity
INCUBATOR/INCUBATOR2	Digital Dry Block Incubators	1
IS050	InSite <i>Salmonella</i> (Environmental <i>Salmonella</i> Species Test)	50
MS2-EB	MicroSnap EB	100
MS2-COLIFORM	MicroSnap Coliform	100
MS2-ECOLI	MicroSnap <i>E. coli</i>	100
MS2-TOTAL	MicroSnap Total	100
See catalog	BAX System Q7	1