

Aller flow Gluten

Hygiena® AllerFlow Gluten is a rapid and convenient test for detecting gluten residue on surfaces as part of an allergen control program. Combining Hygiena's convenient sample collection swab design with classic lateral flow technology, AllerFlow Gluten makes gluten testing easier than ever.

Benefits:

- Convenient and quick sample collection requires minimal training
- Pre-moistened swab enables reliable collection to ensure accurate results
- Simple two-part kit includes all required materials no pipettes, bottles, scales, or heating blocks required
- Pre-measured, self-contained extraction buffer provides consistent results with minimal handling
- Easy-to-interpret lateral flow test results in 10 minutes
- Specific to the toxic fraction of gluten from wheat, barley, rye, and oat
- No cross-reactivity with soy, rice, or corn

Sensitivity:

Sensitive to 5 μg of gluten per swab (5 ppm) within 10 minutes

Positive Result:




Negative Result:




Catalog No.	Item No.	Description	Quantity
ALF-GL	ASY0386	AllerFlow Gluten Kit	25

www.hygiena.com ALG-PS-3004-REVC

Allerflow Gluten Procedure

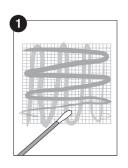

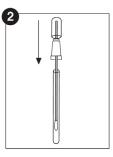









See an instructional demo at HygienaTV or UHygienaTV

Related Products

EnSURE® Touch Monitoring System

EnSURE Touch is the next-generation quality monitoring system. It provides rapid and accurate sanitation verification to support your food safety program. EnSURE Touch features a 5" touch screen, wireless sync technology, and cloud-based software. EnSURE Touch is designed to adapt to your workplace, providing the data you need for complex multi-location supply chain monitoring, audit and risk management, and recall prevention.





AllerSnap® High-Sensitivity Allergen Prevention Test

AllerSnap is a colorimetric allergen prevention test for broad range allergen screening in 15-30 minutes, making it ideal for everyday use.



SuperSnap® High Sensitivity Surface ATP Test

SuperSnap is the most sensitive ATP surface test in Hygiena's ATP hygiene monitoring portfolio of all-in-one devices. Designed to work with all Hygiena luminometers, SuperSnap is an all-in-one test which detects extremely low levels of ATP. SuperSnap is predominantly used in food and beverage processing facilities as an allergen cross contamination prevention.



Catalog No.	Description	Quantity
ETOUCH	EnSURE Touch Monitoring System	1
ALS-100	AllerSnap High-Sensitivity Allergen Prevention Test	100
SUS3000/SUS3000X	SuperSnap High-Sensitivity Surface ATP Test	100

www.hygiena.com ALG-PS-3004-REVC